Skin Cool Pat

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TE800TW

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
Thank you for your purchase.
Before using this product, be
sure to read the manual and
use the product correctly.
After reading, keep the manual
in a location where users of
the product can find it easily.

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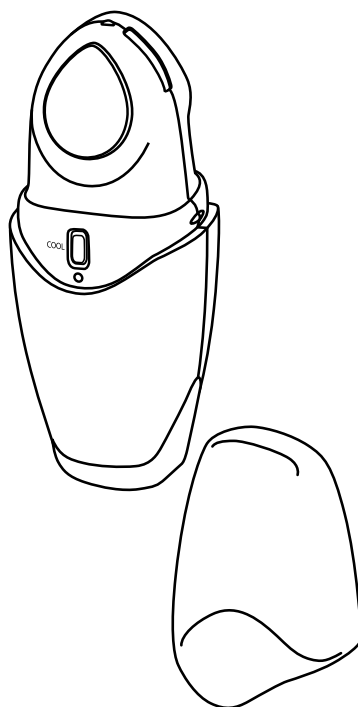Before Using 5

How to Use 6-8

When to Charge/
When Using Overseas/
When Using while Plugged In 8

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『Inspecting Skin Cool Pat that has been used for a long time!』

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●Make sure that the Cool/LED Head parts are undamaged before use.

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Safety Warnings

- Make sure to read these safety warnings before using, and use the product correctly.
- The warnings listed here are to ensure the product is used properly and safely, and avoid any injury to the user or others or damage to property in advance. Make sure to follow the warnings.
- Safety warnings are separated into categories as follows.

Danger

For actions that have a significant chance to cause serious injury or death if done improperly.

Alert

For actions that could cause serious injury or death if done improperly.

Warning

For actions that could cause injury or property damage if done improperly.

Examples of Symbols

This symbol indicates an action that must not be performed.

(The image to the left indicates "no taking the product apart.")

This symbol indicates an action that must be performed.

(Always perform an action indicated by the image to the left)

Alert

Consult with a doctor before use if any of the following apply.

- Heart problems
- Acute (painful) illness ● Fever
- Abnormal blood pressure ● Infectious disease
- Tuberculosis ● Malignant tumor
- Facial neuralgia ● Pregnancy
- While on period ● Atopic eczema
- Allergies ● Particularly sensitive skin
- Hemophilia ● Receiving dental treatment
- Capillary dilatation due to long-term steroidal hormone treatments or liver dysfunction
- While under the influence of drugs or alcohol, or extremely tired

There is a risk of accident or injury.

Required

Do not use on the following areas.

- Areas where plastic surgery was performed
 - Eyeballs and eyelids ● Injuries ● Facial melanosis
 - Inside the mouth or other mucous areas
 - Areas inflamed with pimples or other breakouts
 - Areas with impairment of pain sensation or of perception
 - Areas with metal, plastic, silicone or other objects under the skin
 - Areas that are itchy, have a rash or other stains with physical stimuli
 - Areas inflamed due make-up or other reasons
- There is a risk of accident or injury.

Prohibited

Do not use the product on wounds or inflammation of skin.

It can cause skin problems or pain.

Do not use the product for uses other than those listed in the manual.

It can cause accidents or malfunction, or damage the product.

Prohibited

Do not let children use the product.
Do not use in reach of toddlers.

There is a risk of accident or injury.

Place out of reach of toddlers.

The cap and other small parts may be swallowed accidentally.

Do not use on toddlers, children or pets.

There is a risk of accident or injury.

There is a risk of the power cord disconnecting when used improperly or for long periods, so do not use the product in the following situations:

- If the power cord is twisted and does not return properly.
 - If the base of the power cord or adapter has warped or is damaged.
- If the power cord is disconnected, there is a risk of sparks occurring and causing burns or electric shock.

Do not use when the power cord is damaged or loose in the socket.


There is a risk of fire or electric shock.

Do not damage the power cord. (Do not pull, twist, or bend it unnecessarily, do not place heavy objects on it, do not place it near heat, do not bundle it, do not squeeze it, do not modify it)

There is a risk of fire or electric shock.

Do not wrap the power cord around the power adapter.

There is a risk that the power cord will disconnect and cause fire, electric shocks or burns.

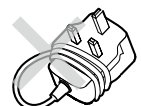Safety Warnings

Prohibited

Do not use with power sources other than AC100-240V.

There is a risk of fire or electric shock.

Do not look at the LED Care light directly for long periods of time.

There is a risk of accident or injury.

If you experience pain or discomfort while using the product or issues with your skin after using the product, stop using the product and consult a medical specialist.

It could cause skin problems.

Use Cool Facial Care with an ambient temperature between 10 and 40°C.

There is a risk of skin damage.

Required

If the product makes abnormal noises or some other irregularity occurs, stop use immediately.

There is a risk of accident or injury.

Follow directions in manual for frequency and duration of use.

There is a risk of accident or injury.

When first using the product, track condition of skin.


There is a risk of skin problems or poor health.

Remove Power Adapter

Remove the power adapter from the outlet when not in use.

There is a risk of fire or electric shock.

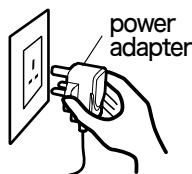No Wet Hands

Do not plug in or remove the power adapter with wet hands.

There is a risk of malfunction, electric shock or short circuit.

No Getting Wet

Do not wash with water. Do not place the product in water or splash it with water.

There is a risk of malfunction or short circuit.

No Wet Areas

Do not use or store the product in a place with high moisture or where it might get wet, such as a bathroom.

There is a risk of malfunction or short circuit.

No Dismanting

Do not modify the product. Do not dismantle or repair the product if you are not a repair technician.

There is a risk of fire, electric shock or injury.

Warning

Do not use at the same time as other beauty goods.

There is a risk of accident or injury.

Do not drop or hit the product.

Do not apply strong impacts to the product.

There is a risk of electric shock, fire or injury.

Do not wipe the product with benzene, paint thinner or alcohol, nor spray disinfectant on the product.

There is a risk of cracking, electric shock or fire.

Prohibited

Do not use power adapters other than the one that comes with the product. Do not charge other products with the power adapter.

There is a risk of electric shock or fire.

Do not use for commercial purposes. Do not put it in a place unknown numbers of people may use it.

Continued use can cause malfunctions.

Do not use if damaged.

There is a risk of electric shock, fire or injury.

Do not share with family or strangers.

There is a risk of infection.

Before using, check that the Cool Head and LED Head are not damaged or deformed.

There is a risk of damaging skin.

Required

Remove contact lenses and accessories such as earrings and necklaces before using.

There is a risk of accident or injury.

When removing the power adapter from the outlet, grasp the adapter itself instead of the power cord.

There is a risk of fire or electric shock.

Names and Functions of Parts

Cool Head

Place on the skin and use for Cool Facial Care.

Cool Switch

Pushing the switch turns it ON and causes the Cool Facial Care power light to blink. Pushing again turns it OFF and turns off the Cool Facial Care power light. Pushing the LED Switch while the Cool Switch is on will switch to LED Care and turn off the Cool Facial Care power light.

Cool Facial Care Power Light

Blinks when the Cool Switch is pushed, and becomes solid when cooling preparation is complete.

Battery Light

[Normal Use]
Off: Charge remaining
Blinking: Low charge
[While Charging]
On: Charging
Off: Charging complete

Air Intake

Do not put objects in here.

LED Head

Place on the skin and use for LED Care.

LED Switch

When the switch is pushed, it turns ON and the LED Head blinking pattern begins. Push again to turn OFF and the LED Head turns off. If the Cool Switch is pushed while in use, the LED Head will turn off and it will switch to Cool Facial Care.

Nickel Metal Hydride Battery (Internal)


Air Exhaust

Do not put objects in here.

Main Body Socket

■Accessories

Power Adapter

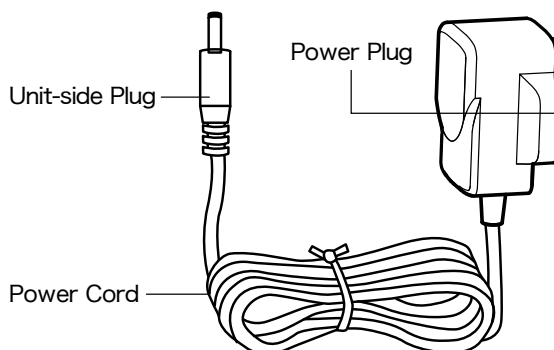Cap

*Follow regional laws when disposing of this product.

How to Charge

1 Confirm the Cool Switch and LED Switch are OFF, and connect the power adapter to the main unit and outlet.



- Pushing the unit-side plug of the power adapter too hard or forcefully bending it can damage the unit socket.

2 When the Battery Light turns off, unplug the power adapter from the main unit and outlet.

- The average charging time is approximately 5 hours.
- While charging, the main unit and power adapter may become warm, but this is not a malfunction.
- One charge allows for use of approximately 5 sets. (1 set: 5 minutes of LED Care and 5 minutes of Cool Facial Care)

(Based on factors such as conditions of use and maintenance, the number of sets may decrease)

The Battery Light will turn on while charging, and turn off when charging is complete.

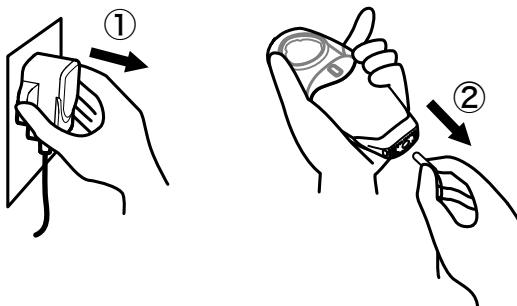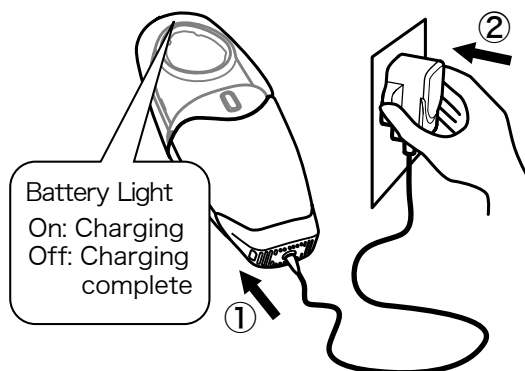How to Have a Long-lasting Nickel Metal Hydride Battery

- When storing the product, fully charge it before storage.
- When not using for a long time, charge periodically every six months.

When the operating time after charging becomes extremely short, the following procedure can refresh the battery and improve operating time.

1. Push both the Cool Switch and LED Switch simultaneously for over 10 seconds. When the Cool Facial Care power light and Battery Light start blinking, stop pushing the switches. Let the charge run out until the lights turn off.
2. Charge as normal.

Warning

- Once charging is complete, remove the power adapter from the outlet.
- Use with an ambient temperature between 10 and 40°C.

Before Using

Alert

Consult with a doctor before use if any of the following apply.

There is a risk of accident or injury.

- Heart problems
- Acute (painful) illness
- Fever
- Abnormal blood pressure
- Infectious disease
- Tuberculosis
- Malignant tumor
- Facial neuralgia
- Pregnancy
- While on period
- Atopic eczema
- Allergies
- Particularly sensitive skin
- Hemophilia
- Receiving dental treatment
- Capillary dilatation due to long-term steroidal hormone treatments or liver dysfunction
- While under the influence of drugs or alcohol, or extremely tired

Do not use on the following areas.

There is a risk of accident or injury.

- Areas where plastic surgery was performed
- Eyeballs and eyelids
- Injuries
- Facial melanosis
- Inside the mouth or other mucous areas
- Areas inflamed with pimples or other breakouts
- Areas with impairment of pain sensation or of perception
- Areas with metal, plastic, silicone or other objects under the skin
- Areas that are itchy, have a rash or other stains
- Areas inflamed due make-up or other reasons

If you experience any abnormalities with your skin or body while using, immediately stop use and consult with a medical specialist.
There is a risk of accident or poor health.

*** Guidelines for skin care are to use each type of care twice per day, for no more than 5 minutes each time.**

Using continuously for over 5 minutes can lead to skin problems, so please refrain from doing so.

This appliance is not designed to be used by people (including children) with reduced physical, sensorial or mental capacities, or by people with no prior experience or knowledge, except where they have been able to benefit from supervision or previous instruction relating to use of the appliance from a person responsible for their safety.

Wash your face and clean off any dirt.

- Use cleanser to remove makeup.

Use on skin without any makeup.

For dry skin, apply beauty wash before use.

After applying beauty wash, be sure skin is dry before using.

*Do not use milky or gel beauty washes.

How to Use

This product has 2 kinds of care functions.

LED Care: Using light to improve skin elasticity

Cool Facial Care: Cooling to tighten pores

Do not use either Care more than twice per day (5 minutes or less each time).
The LED Switch and Cool Switch automatically turn off 5 minutes after being pushed.

Step 1 LED Care

Alert

- When using around the eyes, close your eyes as it is bright.
 - Avoid looking directly at the light for long periods of time.
- There is a risk of accident or injury.


Use for 5 minutes or less.

1 Remove cap.

2 Push down on the LED Switch for a while.

- The LED Head will start to slowly blink.

LED blinking
pattern
(Image)

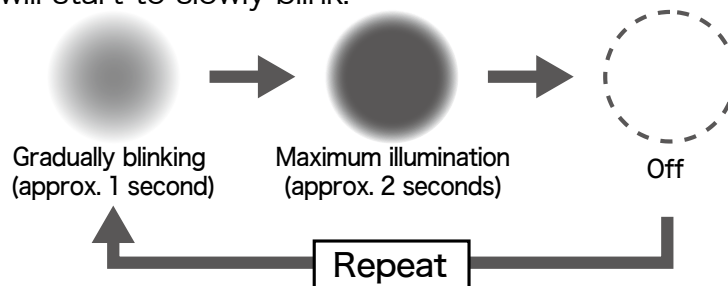3 Gently apply the Head mainly to the area where you have issues with your pores (nose, cheeks, etc.) and slowly move it.


- Use for 5 minutes or less. The power will automatically shut off 5 minutes after pushing the LED Switch.

4 After 5 minutes the power will automatically shut off.

- When that happens, a beeping noise will alert you.

When you want to stop partway through

- Push the LED Switch again and the LED Head will turn off.

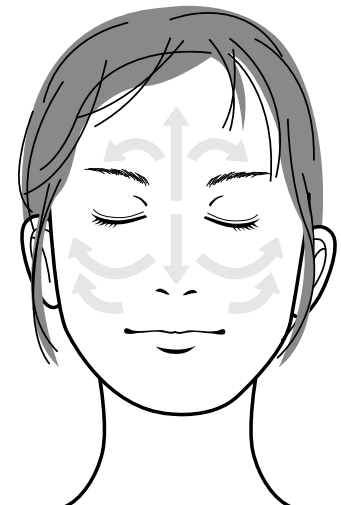Step 2 Cool Facial Care

How to Use

⚠ Alert	Do not use near eyes or ears. (There is a risk of accident or poor health)
	Do not use when skin is injured or cut. (There is a risk of skin problems or pain)
⚠ Warning	Use for 1 second or less in the same area. (There is a risk of skin problems or pain)
Use for 5 minutes or less (including cooling preparation time).	

1 Remove cap.

2 Push down on the Cool Switch for a while.

- The Cool Facial Care power light will start to blink.
- Wait until the cooling preparation ends and the Cool Facial Care power light becomes steady. (Preparation takes approximately 60 seconds)

3 Begin using once the Cool Facial Care power light becomes steady. Check the Cool Head with your fingertip.

- If the Cool Head feels too cold, use it in an area with a high ambient temperature.

The noise while using it is the cooler. It is not a malfunction.

4 Hold so as to not block the Air Intake or Exhaust.

- If the Air Intake or Exhaust are blocked, the Head will not cool properly.

5 Gently apply the Head mainly to the area where you have issues with your pores (nose, cheeks, etc.) and slowly move it while patting.

- For use on the same area, apply the head for 1 second and remove it for 1 second, repeatedly.
- The power will automatically shut off 5 minutes after pushing the Cool Switch.

The main unit will warm up while in use, but this is not a malfunction.




6 After 5 minutes the power will automatically shut off.

When that happens, a beeping noise will alert you.


When you want to stop partway through

- Push the Cool Switch again and the Cool Facial Care power light will turn off and cooling operation will end.

After turning off the power, the temperature of the Cool Head and main unit will temporarily increase, but this is not a malfunction.

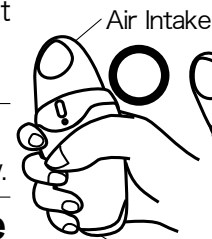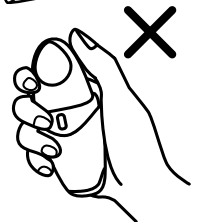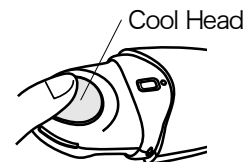Air Exhaust

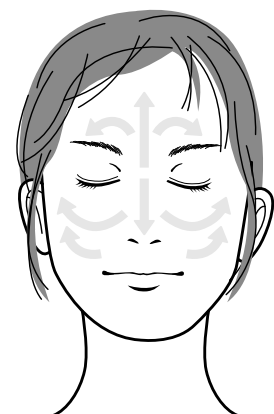Warning

- Do not leave in direct sunlight (inside a car in summer, etc.) or near heating devices. This can cause deformation or discoloration.

How to Use

How to Store Long-term

Perform maintenance (Page 9) before storing.

1. Put the cap on the head.
2. Make sure to charge the product.
3. Store it in an area with limited moisture or dust.

【When to Charge】

- The signal for low batteries is when the Battery Light begins blinking while the power is on and a warning beeping occurs for 10 seconds before the power automatically turns off.

→Charge immediately.

To continue using, use it while plugged in.

Battery Light

【When Using Overseas】

- This product can be used with AC100-240V power.

【When Using while Plugged In】

1 Connect the power adapter to the main unit and power outlet.

2 Being LED Care or Cool Facial Care.

Follow directions (Pages 6-7) for use.

3 After use, remove the power adapter from the outlet and main unit.

To continue charging,



follow the directions on How to Charge (Page 4).

Care and Maintenance

Make sure the power is OFF and the power adapter is not attached before performing maintenance.

Use a soft, dry cloth to wipe off dirt on the main unit.

- When dirt does not come off easily
Soak a dry cloth in soapy water or a neutral cleaning agent cut with water, and wring it out thoroughly before wiping with it.
- Periodic maintenance
Periodically (at least once a month) clean the Air Intake and Intake and grooves around the heads.
Use a soft toothbrush to remove dirt and dust.

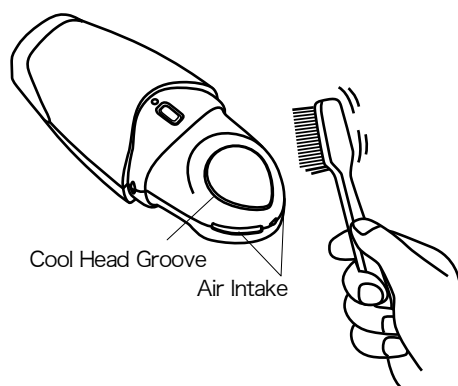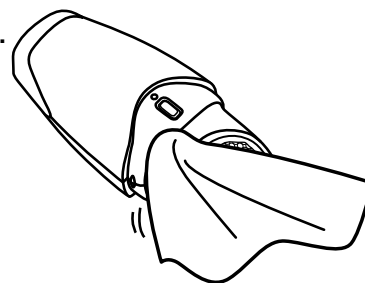*Make sure dust or dirt does not enter the main unit through the Air Intake.

Warning

- Make sure the power is OFF. If you accidentally turn on the power, it could shock you or cause a malfunction or short circuit.
- Do not wash the main unit with water. Do not immerse it in water or splash water on it. It could cause malfunctions.
- Do not wipe with benzene, paint thinner or alcohol. It could cause deformation or discoloration.

Specifications	Product Name		Skin Cool Pat	Rated Time	LED Care: 5 minutes (auto OFF) Cool Facial Care: 5 minutes (auto OFF)
	Product Number		TE800TW		
	Power Source	Power Adapter	AC100-240V 50/60Hz	Charge Time	Approx. 5 hours
		Main Unit	Nickel Metal Hydride Battery (DC3.6V, 750mAh) 1 unit	Weight	120 g (main unit)
	Power Consumption	Cool Facial Care	When plugged in 10W		
		LED Care	When plugged in 2W	Accessories	Cap Power Adapter
	Dimensions		Height 111 mm x width 46 mm x depth 40 mm (main unit)		

Troubleshooting

Issues	Possible Reasons	Solution
The Cool Head is not cooling.	The Cool Switch is not ON.	Push the Cool Switch and confirm the Cool Facial Care power light changes from blinking to a steady light.
	The Air Intake or Exhaust is blocked.	Hold properly to ensure the Air Intake and Exhaust are not blocked.
	There is no charge left.	Charge it or use it while plugged in.
The LED Head isn't lighting up.	The LED Switch is not ON.	Turn on the LED Switch.
	There is no charge left.	Charge it or use it while plugged in.
The Battery Light starts blinking too soon.	It was not used for over 3 months. (Natural discharge)	Charge it.
It can only be used a few times after charging.	The product is charged despite there still being charge remaining. (Top-up charge)	Use the battery fully before charging.
	The nickel metal hydride battery is getting old.	Contact the information number on the product box.
Charging does not work.	The power adapter is not connected properly.	Make sure the power adapter is securely connected to the main unit and the outlet.
	The power adapter is malfunctioning.	Contact the information number on the product box.
The power cord or power adapter are abnormally hot.	The power adapter is not connected properly.	Firmly connect the unit-side plug of the power adapter to the main unit socket.
	The plug is not connected to the outlet firmly.	The issue may be with the outlet. Contact the information number on the product box.
It cannot be used with the power adapter plugged in. (Not turning on while plugged in)	The power adapter is malfunctioning. Or the main unit is malfunctioning.	Contact the information number on the product box.

In the case of other errors or malfunctions, immediately stop use and contact the information number on the product box.

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